



DMS 12-lead Resting ECG Work Station



Blue Tooth Sample Device



Wire Sample Device

Whole Set Include:

Item	Quantity
Blue-tooth kit	1
ECG CABLE	1
SUCKS	1
CLIPS	1
ECB CABLE ADAPTER	1
POUCH	1
Bluetooth ECG Transmitter and Receiver	1
SECURITY KEY	1
Product warranty card	1
CERTIFICATION LABEL	1



Wire Kit 12 lead Resting ECG

Sample device with Premier Software

Sample device with Standard Software


Whole Set Include:

Item	Quantity
ECG Acquisition Box	1
ECG CABLE	1
SUCKS	1
CLIPS	1
ECB CABLE ADAPTER	1
DMS300-BTT02 MANUAL	1
SECURITY KEY	1
Product warranty card	1
CERTIFICATION LABEL	1

Features of Software

Software (Standard Version / Premier Version)


New Computerized Resting ECG 12-Lead ECG with Interval/Amplitudes/Axis measurements/with or without automatic interpretation developed and made by DM Software, implementing all the advanced holter ECG technology from DM Software, and many unique analysis functions:

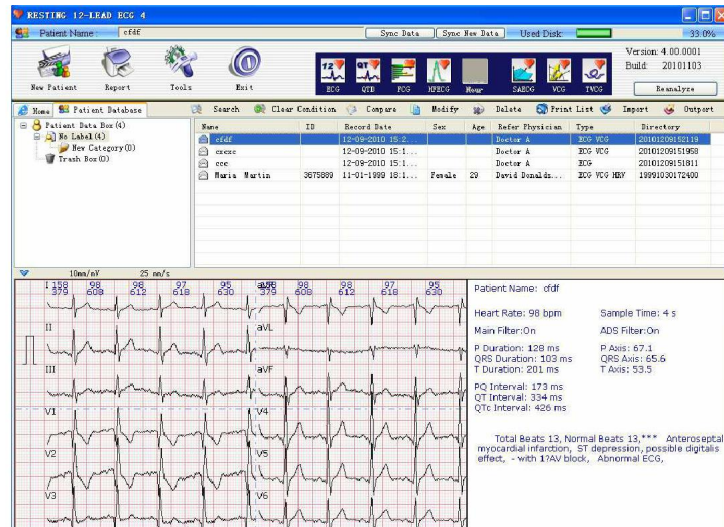


Version 4.0 Standard	12 Lead Resting ECG with or without auto interpretation
	Vector ECG
Version 4.0 Premier	12 Lead Resting ECG with or without auto interpretation
	Vector ECG
	Timed Vector ECG
	SAECG (Signal Averaged ECG) Ventricular Late Potentials analysis
	QTD Analysis
	Hour ECG
	High Frequency ECG (HFECG)
	Frequency Cardiogram (FCG)
Simple Exercise ECG	


Function of Software

1. Interface

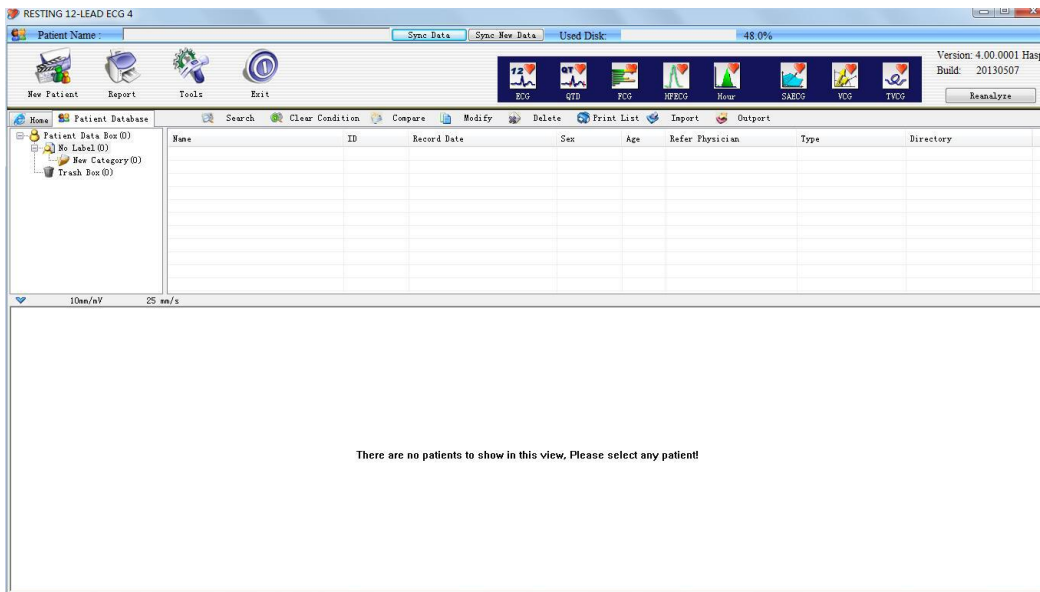


2. New Patient or Demo

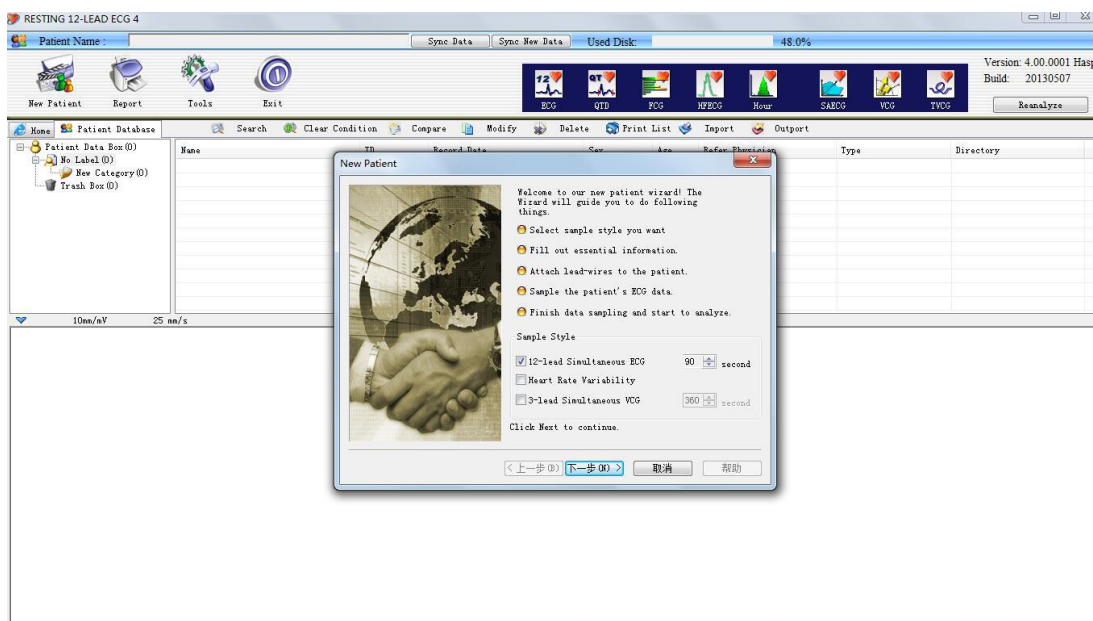

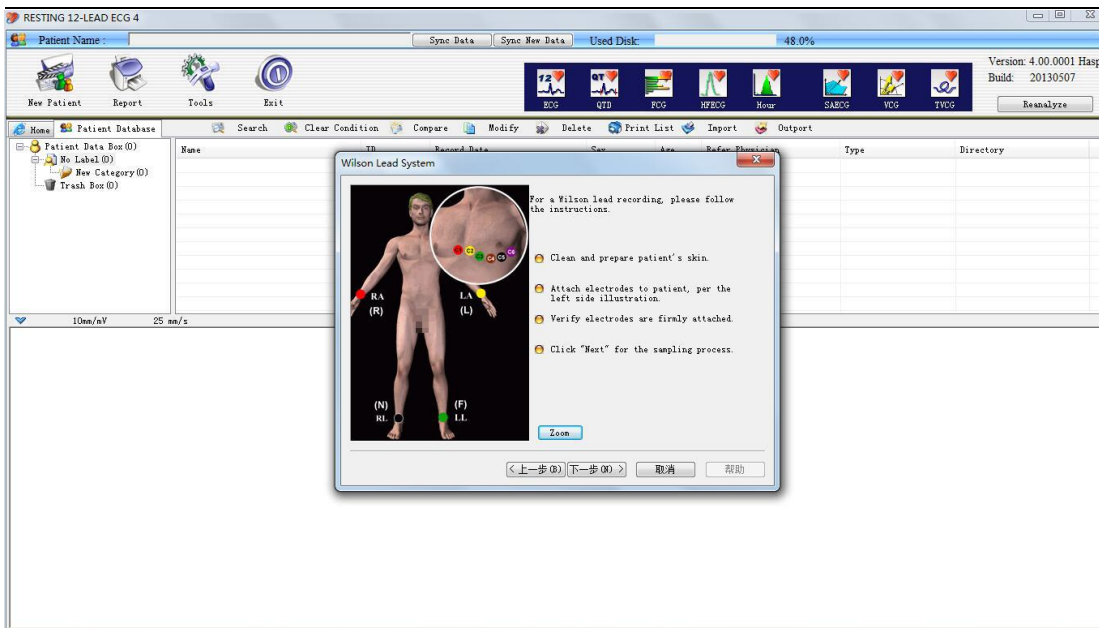


3.Operation and Real Time Monitoring

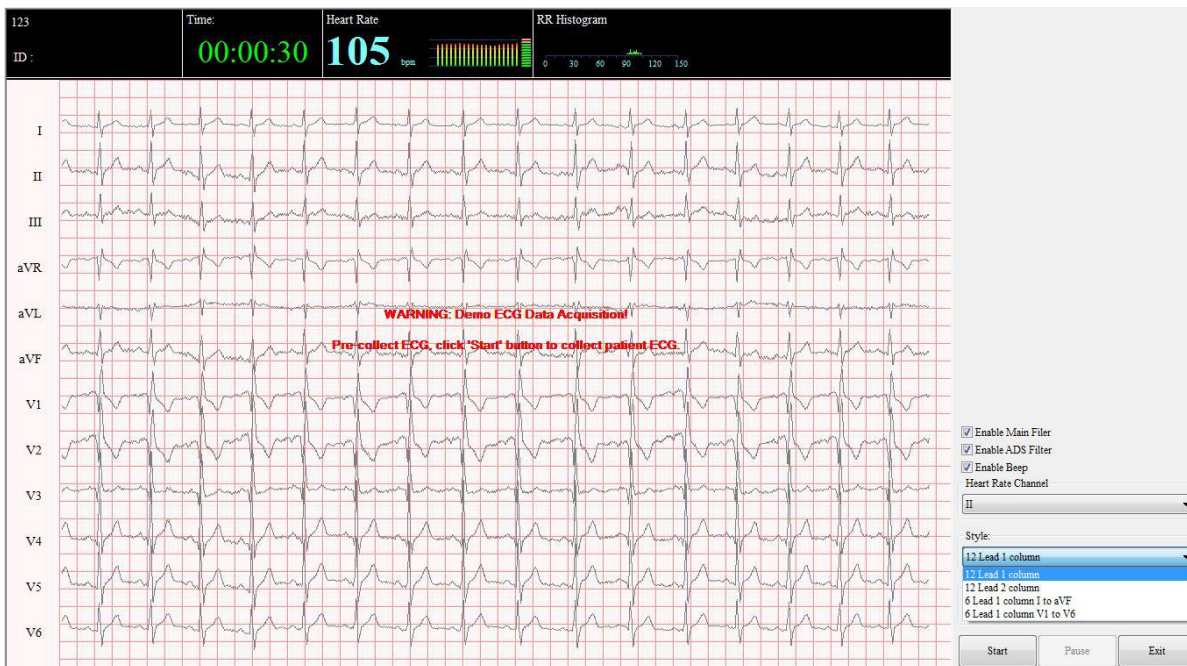


4.Analysis Functions

ECG Review

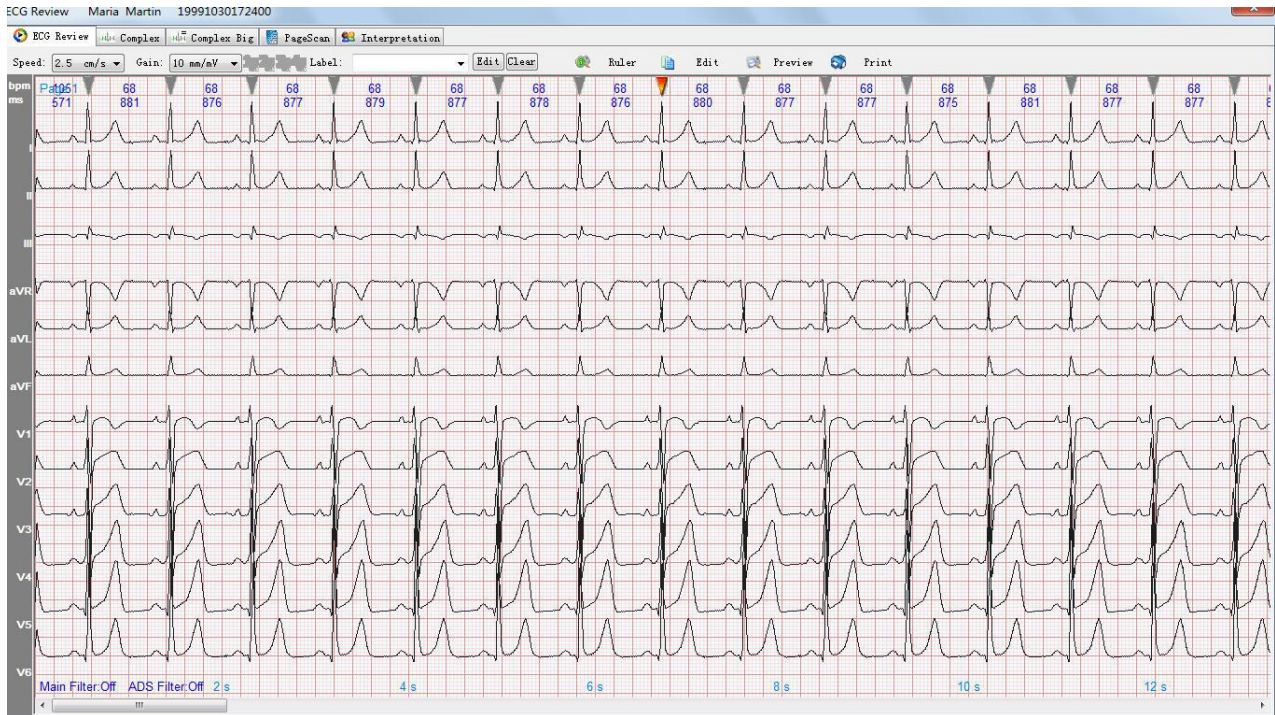

FCG

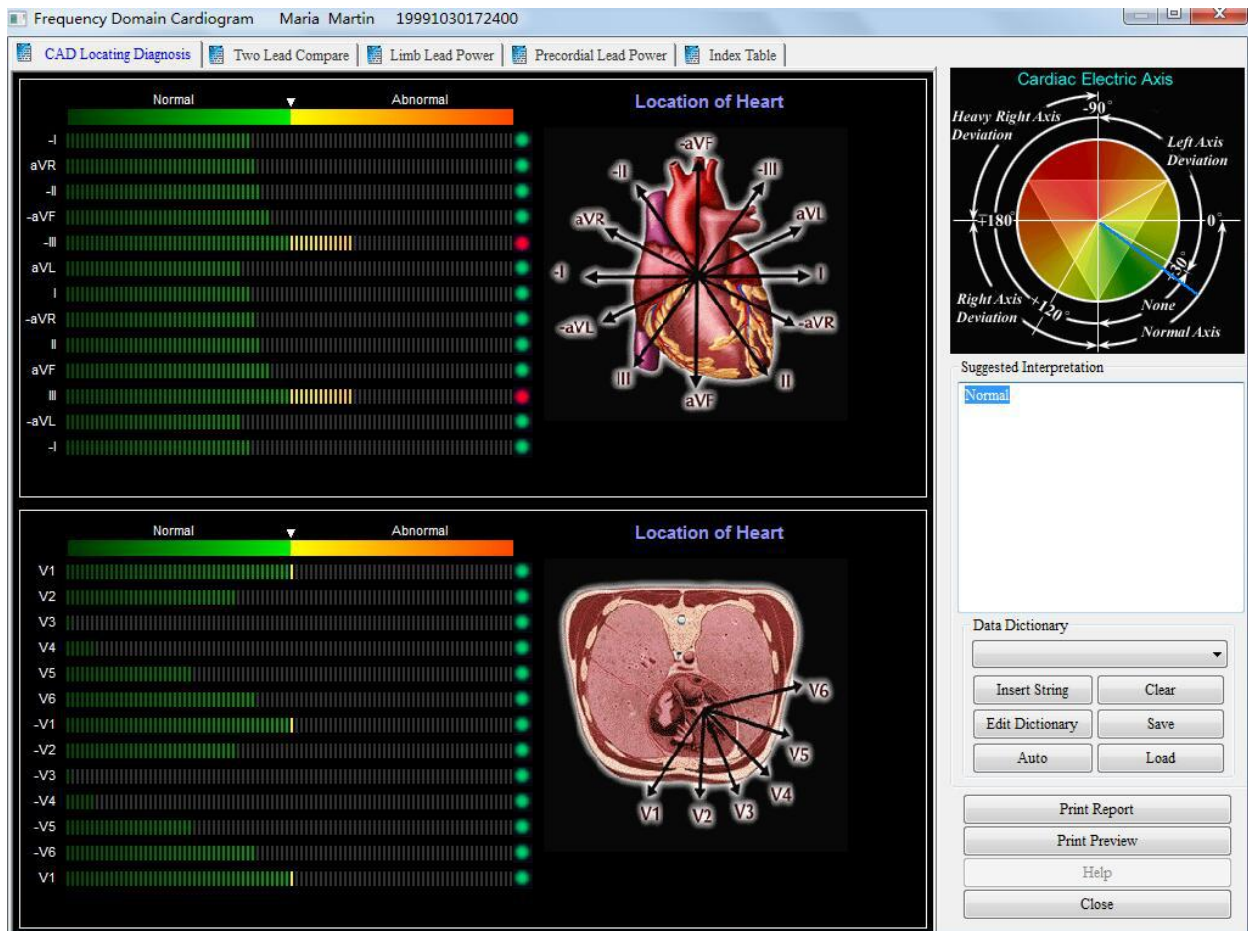





HFECG




QTD

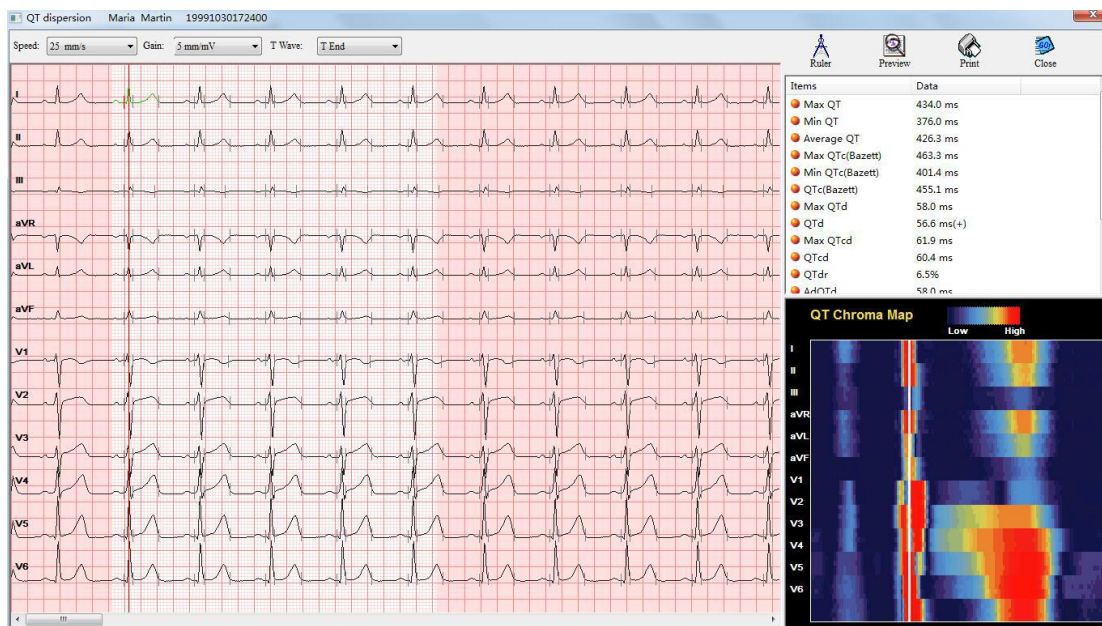


SAECG

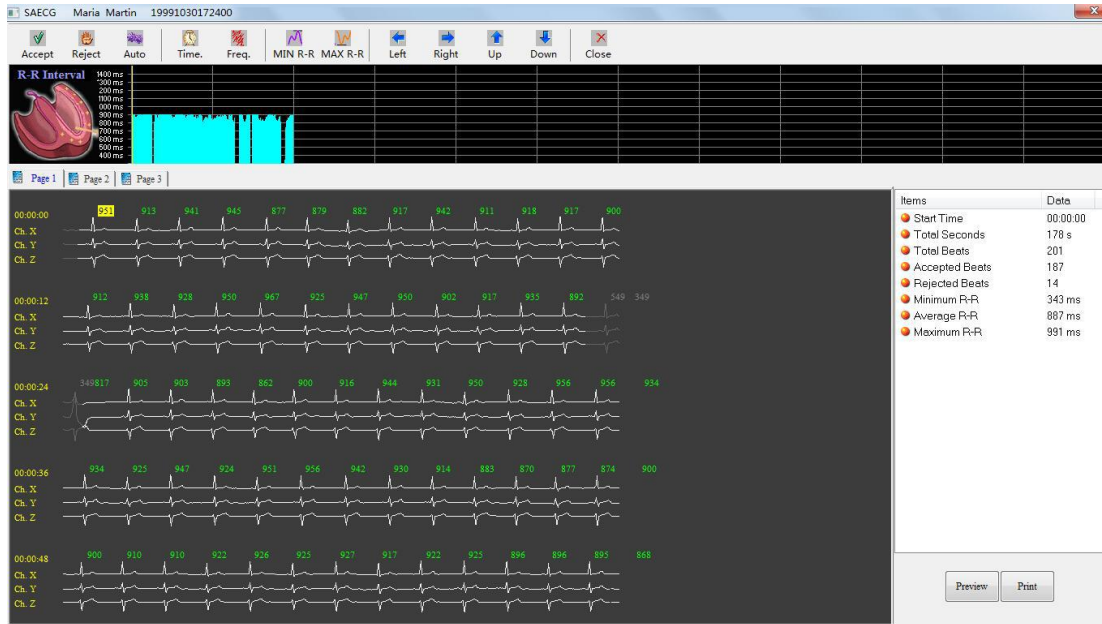


TVCG

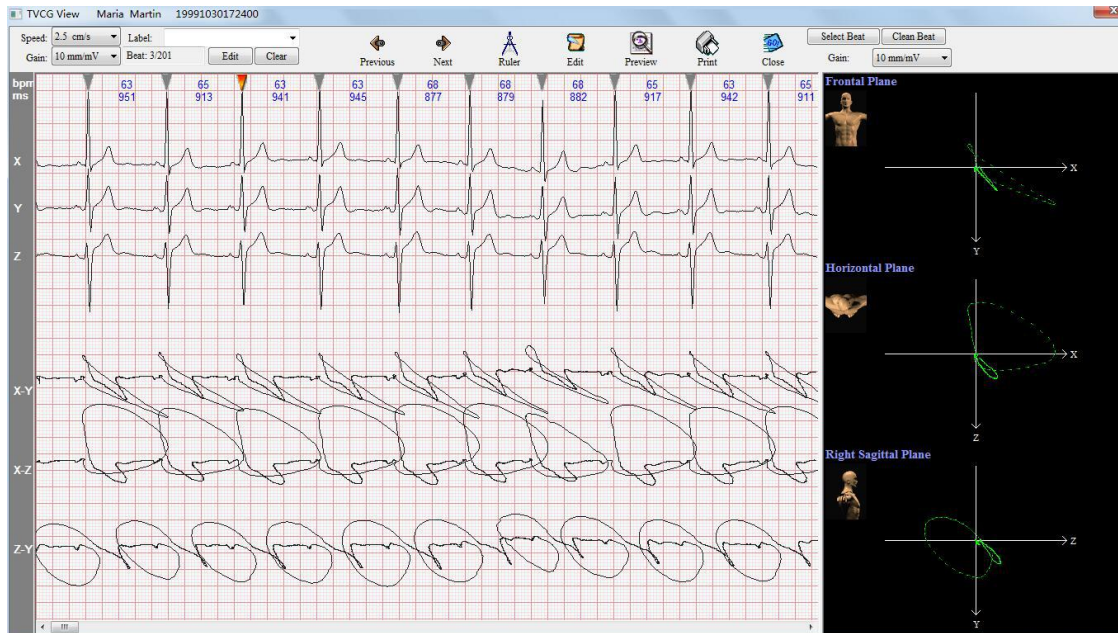

VECG

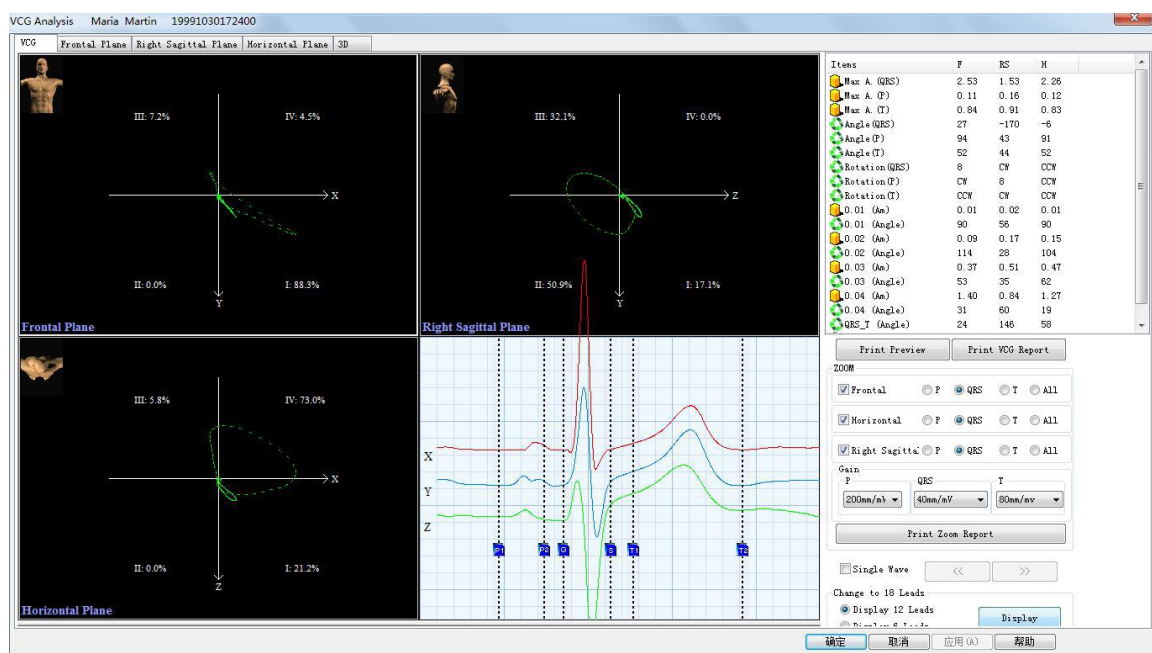

5.Reports

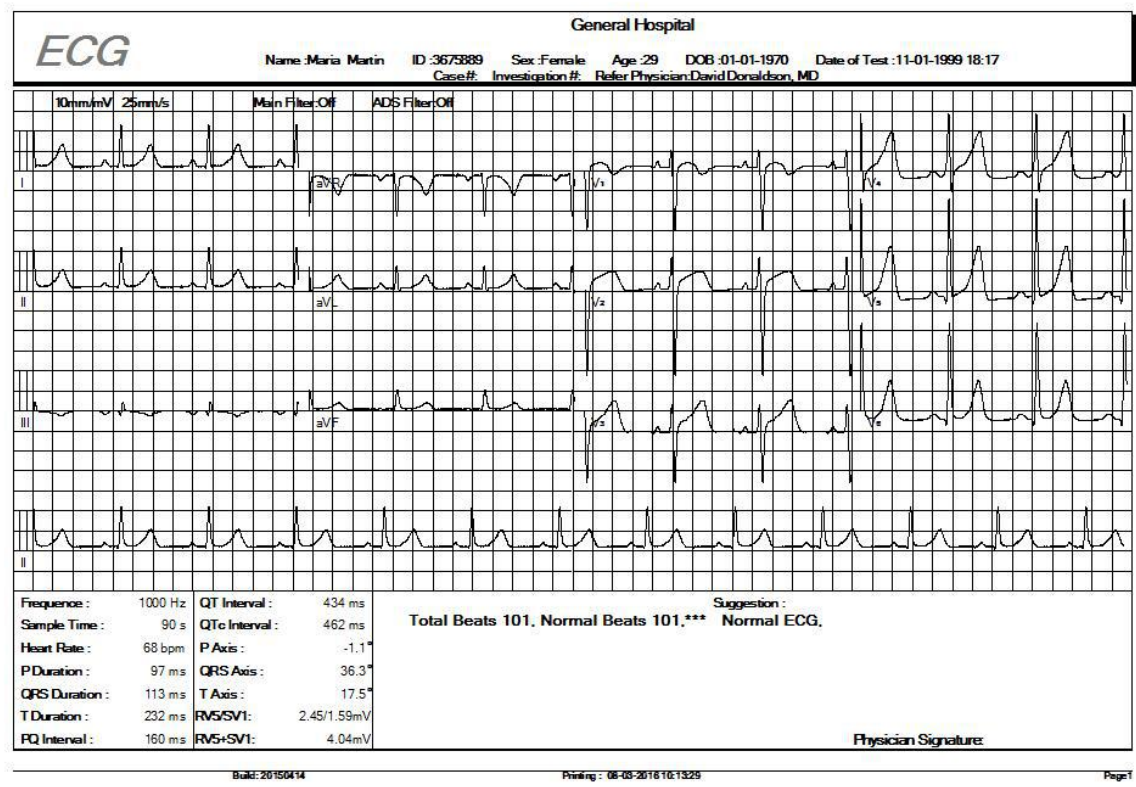

Print Other Reports



Report Maria Martin 19991030172400

ECG Option
 Print 12 Leads ECG
 Suggested Interpretation:
 Total Beats 101, Normal Beats 101,*** Normal ECG,

ECG Report
 Landscape (3*4+2)
 Landscape (12)
 Landscape (3*4+1) with Table
 Landscape (6*2+1) with Tabl.
 Landscape (6*2+1)
 Table Page
 Portrait (12) with Data
 Avg ECG

Other
 Print QTd
 Print FCG
 Print HFECG
 Print Hour ECG
 Print SAECG
 Print VCG
 Print TVCG

Print Selected Send to EMR GDT (0)
 Save ASCII ECG Cancel

QRS Cardiac Electric Axis
Cardiac Electric Axis
 Heavy Right Axis Deviation Left Axis Deviation
 Right Axis Deviation None Normal Axis
 QRS Axis: 38.3 degree

Hardware Description


Bluetooth Kit:





Leads	12-Lead, 10 electrodes
Sampling Rate	12 Bits, 4096Hz
Dimension	Trans. 13.5 x 7.5 x 2.4 cm Receiver. 13.9 x 12 x 3 cm
Weight	Trans. 110 g Receiver 240 g
SAECG * & VCG *	Yes
FCG *	Yes
Power	Trans. 4 x AA battery. (6.0V) Receiver. USB Port (5.0V)
ECG Download	Blue Tooth Transmission 2,4 GHz 100 M
Compatibility	CardioScan Resting 12-Lead ,CardioScan Stress Test

Wire ECG Kit:

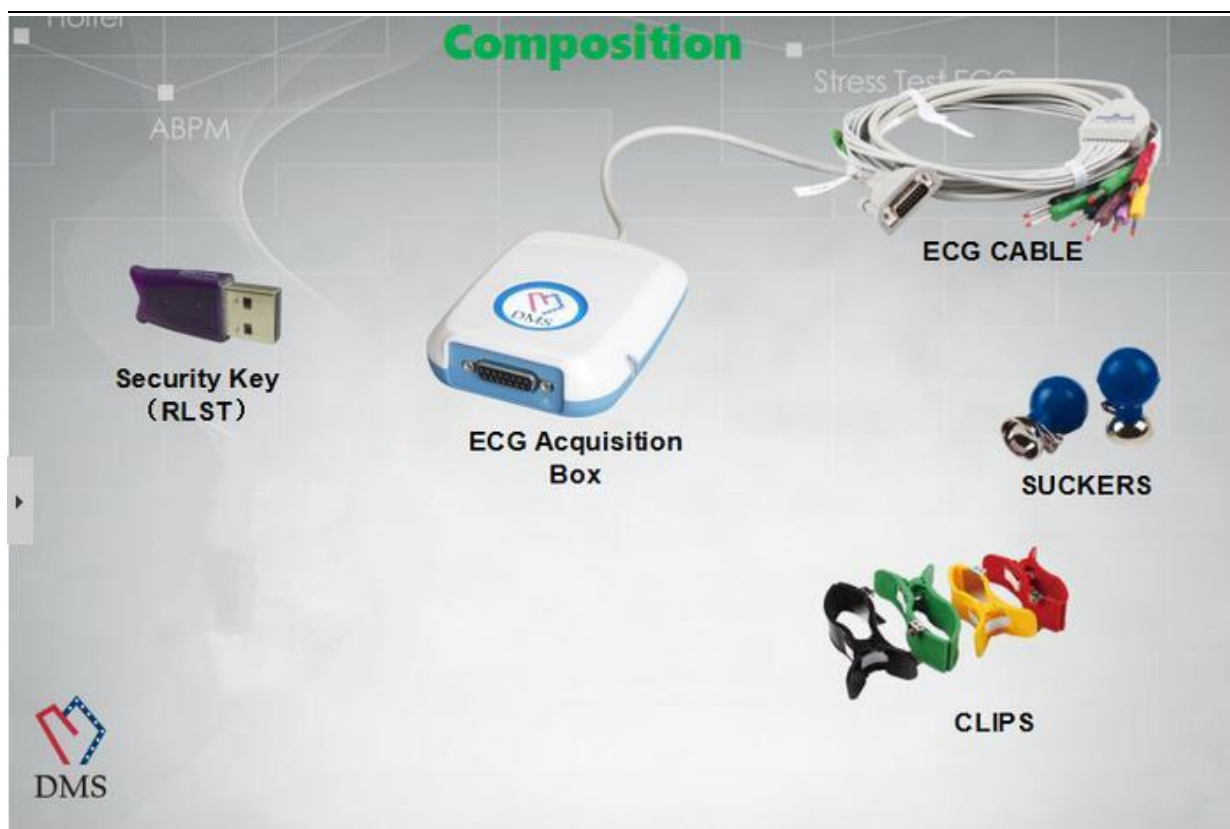

Leads	12-Lead, 10 electrodes
Sampling Rate	16 Bits, 10240Hz
Dimension	12*8.5*2.8 cm
Weight	190g
SAECG * & VCG *	Yes
FCG *	Yes
Power	USB Power Supply,DC 5.0V(+5%,-10%)
Compatibility	CardioScan Resting 12-Lead,CardioScan Stress Test